




PHOTOELECTRIC COLORIMETER

AP-1101 **MADE IN JAPAN**

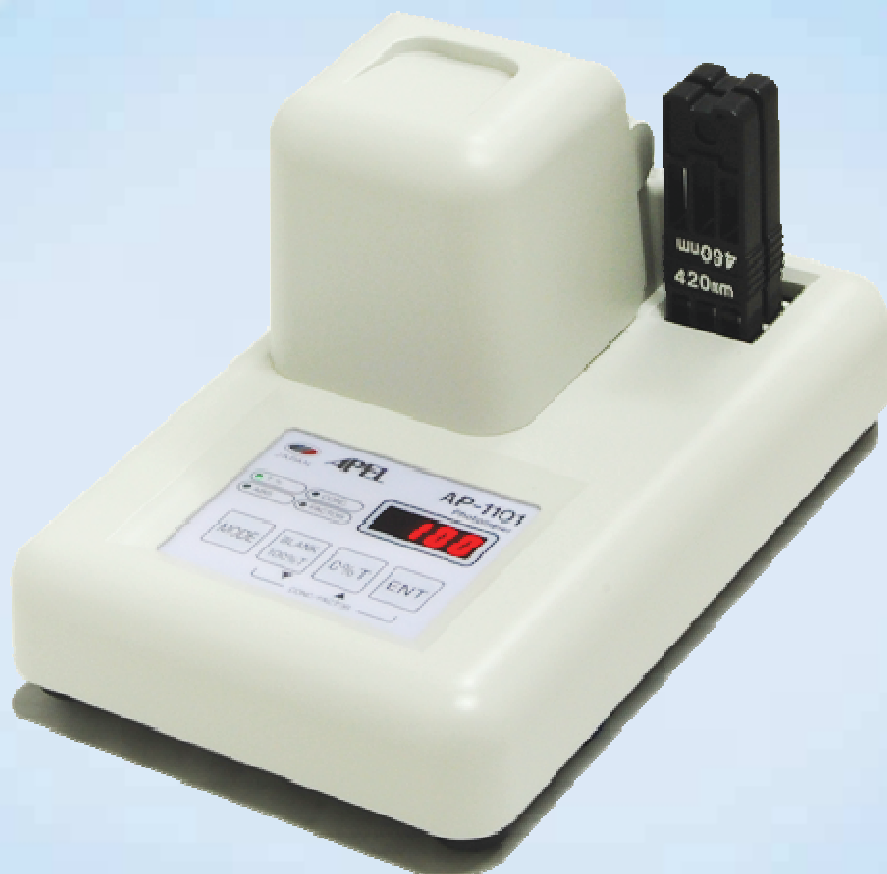

Easy and simple to use

■ BRIGHT INTENSITY LED

High efficiency and bright intensity LED is adopted as the light source, therefore it is sensitive and is not necessary to be changed.

■ DIRECT READING OF CONCENTRATION

By setting the concentration with standard value, the concentration value of each sample can be read directly.

■ SQUARE AND ROUND CUVETTES USABLE

In addition to square cuvettes, round cuvettes can also be used only by fitting out with a supplied adapter.

■ THE SMALLER SAMPLE VOLUME

The minimum sample volume is only 1.0ml.



PHOTOELECTRIC COLORIMETER AP-1101



SPECIFICATIONS:

Display	3 1/2 digit LED				
Built in Filters	420	460	510	540	600
Photometric range	0 - 10%T 0 - 1.999A 0 - 1999C (0-1999F)				
Sample volume	1.0ml (minimum)				
Light source	High bright LED				
Auto zero	Yes				
Power	AC adapter				
Dimensions	150mm(W) x 230mm(L) x 135mm(H)				
Weight	Net 0.65kg (Main body)				

STANDARD ACCESSORIES:

Description		Quantity
Main body		1unit
Filter(Three kinds)	420nm,460nm	1pc
	510nm,540nm	1pc
	600nm	1pc
Measurement cuvette	Square cuvette	2pcs
	Round cuvette	2pcs
Vinyl dust cover		1pc
AC/DC power code		1pc
Operation		1copy

• Specifications and/or appearance are subject to change without any prior notice.



APEL Co.,LTD.

URL: <http://www.apel.co.jp/>