



APEL CO., LTD. JAPAN

**APEL Photoelectric colorimeter
MADE IN JAPAN**

AP-120

Compact Photoelectric Colorimeter

**State-of-the-Art
Colorimeter**



Compact and simple

Built-in Filter

Five filters are built-in.

Auto-Zero by microprocessor control

T=0% and T=100% are automatically set with just a touch of key.

Compact & light weight

AP-120's foot print is small and weights only 0.7kg.

No lamp replacement

AP-120 does not need a lamp replacement thanks to the LED Light source.


**Stray light
Free**



Cover prevents light penetration from outside.

AP-120

Compact Photoelectric Colorimeter

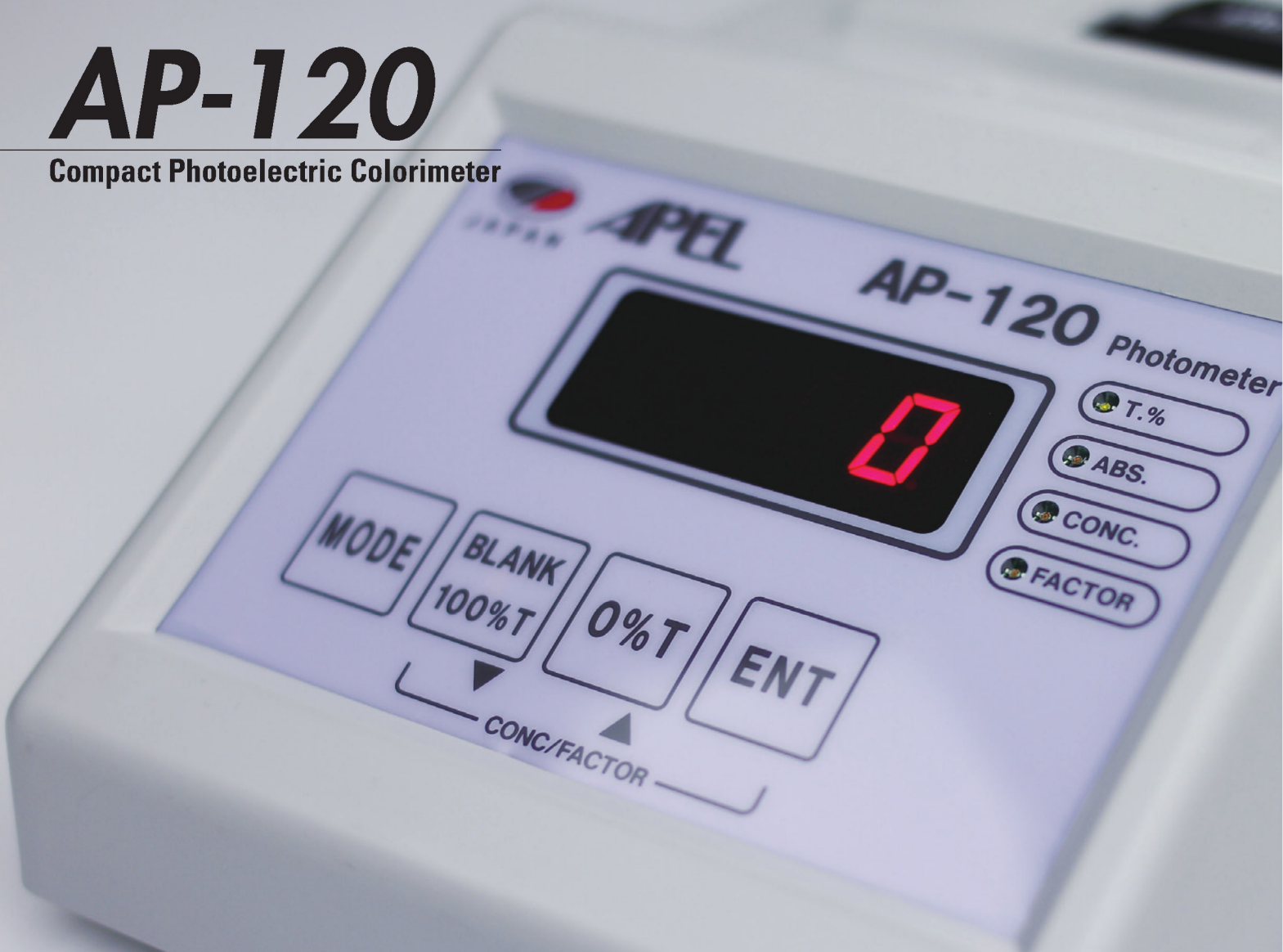

Specifications:

Display	3 1/2 digit LED				
Built in Filters	420	460	510	540	600
Photometric range	0 – 100%T				
	0 – 1.999A				
	0 – 1999C (0 – 1999F)				
Sample volume	1.0 ml (minimum)				
Light source	Bright intensity LED				
Auto zero	Yes				
Power	AC adapter(DC 12V)				
Dimensions	117mm(W)x185mm(L)x77mm(H)				
Weight	Net 0.7kg (Main body)				

Specifications and appearance are subject to change without notice.

Standard Accessories:

Instruction manual
AC adapter(Apel's version)
10×10×45mm square cuvette × 2

MADE IN JAPAN



APEL CO., LTD. JAPAN

Website

<http://www.apel.co.jp>

Distributor